

OTOLOGY PRODUCTS





OTOCELL® EAR WICKS

Description

OTOCELL® Ear Wicks are constructed of a unique medical grade sponge material. The wicks are supplied fully compressed for ease of insertion. The material will expand as fluid is introduced, giving the surgeon time for accurate positioning. OTOCELL® Ear Wicks are available in 8 mm diameter and two lengths (15 mm and 24 mm), with or without fenestration.

Indications for use

For use in the treatment of otitis externa.

Supplied

Sterile

50 wicks per dispenser carton

(1 wick per pouch, 50 pouches per carton)



FEATURES

- Lint and fiber-free sponge material
- Superior liquid absorption and wicking
- Supplied compressed for ease of insertion
- Extremely soft and pliable cell structure

Product Codes

SIZE	REF	
	not fenestrated	fenestrated
15 mm	EW-0815	EW-0815F
24 mm	EW-0824	EW-0824F



OTOCELL® Ear Wick, not fenestrated



OTOCELL® Ear Wick, fenestrated



NEW
Improved Quality

More sponge products are available.

Contact us for your copy of our catalog with RHINOCELL® PVA sponges for Rhinology.



OTOCCELL® INSTRUMENT WIPES

Description

OTOCCELL® Instrument Wipes are 3.25" x 3.25" sheets made of a unique medical grade sponge material. The lint and fiber free material is used for intraoperative removal of debris from delicate optics or instruments which will not snag on its smooth surface. Microsurgical instruments used in tympanoplasty and stapedioplasty can be cleaned safely and efficiently without a risk of introducing any fiber into the sensitive middle ear area.

Indications for use

For intraoperative cleaning of delicate instruments.

Supplied

Sterile

20 wipes per dispenser carton

(1 wipe per pouch)

Dimensions & Product Codes

REF	Dimensions (cm)		
	Thickness	Length	Width
WP-881	0.1	8.25	8.25
WP-882	0.2	8.25	8.25



FEATURES

- Lint and fiber-free sponge material
- Superior liquid absorption and wicking
- Supplied compressed for ease of insertion
- Extremely soft and pliable cell structure

New ✓

OTOLOGICAL SILICONE STRIPS

Description

Silicone Strips, made of medical grade silicone, are used in the process of packing the external auditory canal following surgical procedures such as tympanoplasty or stapedectomy. The strips are intended for temporary use (up to 30 days) and have rounded edges to minimize the risk of trauma to the auditory canal. They can be moistened prior to use using sterile saline solution. The blue color of the strips aids in visualization during removal.

Indications for use

For use as a temporary dressing in the external auditory canal following surgery.

Supplied

Sterile

30 strips per carton

(3 strips / pouch, 10 pouches / carton)

Dimensions & Product Codes

REF	Dimensions (mm)		
	Thickness	Length	Width
SS-4506B	0.12	45	6



FEATURES

- Thin, soft and pliable material
- Economical pre-cut strips
- Non-adherent material
- Transparent blue color for increased visualization

OTOLOGY PRODUCTS

OTOSILK™ OTOLOGICAL STRIPS

Description

OTOSILK™ Otolological Strips are used in the process of packing the external auditory canal following surgical procedures such as tympanoplasty or stapedectomy. The strips can be moistened prior to use using either sterile saline or medication solution. The blue color of the strips aids in visualization during removal. The ends of the OTOSILK™ strips have been heat-sealed to prevent unraveling or shedding.

Indications for use

For use as a temporary dressing in the external auditory canal following surgery.

Supplied

Sterile

200 strips per carton

(10 strips / pouch, 20 pouches / dispenser carton)

Dimensions & Product Codes

REF	Dimensions (mm)	
	Length	Width
OT-17510	45	6.5
OT-16610	45	4



FEATURES

- Convenient, economical pre-cut strips
- Non-adherent material
- Blue color for increased visualization
- Heat-sealed edges

OTOSILK™ GRAFT DRESSING

Description

The OTOSILK™ Graft Dressing is a 10 x 15 cm sheet of blue, nonadherent material for use as a carrier for split thickness skin grafts placed in the external auditory canal. The graft is placed on the dressing and trimmed to the desired size and shape. The dressing and graft combination is then placed in the canal. Following the healing process, the graft dressing is removed, leaving the newly healed graft in position.

Indications for use

For use as a temporary dressing in the external auditory canal following surgery.

Supplied

Sterile

10 dressings per carton

(1 dressing/sterile pouch; 10 pouches/dispenser carton)



FEATURES

- 10 x 15 cm sheet, easily trimmed
- Non-adherent material
- Blue color for increased visualization

Product Code

REF	
OT-4601	OTOSILK™ Graft Dressing



DOYLE EAR DRESSING

Description

The patented Doyle Ear Dressing is a postoperative ear dressing used to protect the external ear following surgery. The unique anatomical shape provides for a better fit and increased comfort. The plastic ear cup protects the external ear area while small openings provide ventilation. The soft fabric strap is color-coordinated and fully adjustable.

Doyle Ear Dressings are available for both adult and pediatric patients. Adult models are available in natural and blue colors. Pediatric models are available in natural, blue, green, orange, and pink colors. A pediatric assorted pack is also available containing one green, one orange, two blue, and two pink dressings.

Each Doyle Ear Dressing contains an ear cup with integral foam pad, color-coordinated adjustable strap, anatomically shaped non-adherent dressing, gauze fluff, adhesive pads, and cotton balls. Pediatric dressings include fun stickers that can be applied to the ear cup for decoration.

Indications for use

For use as an external ear dressing following internal and/or external ear surgery.

Supplied

6 dressings per carton
(1 ear dressing per sterile pouch, 6 pouches per dispenser carton)




FEATURES

- Patented anatomical design
- Comfortable, adjustable strap
- All accessories included
- 4 fun colors and stickers for children
- Adult models in natural or blue

Product Codes

REF	
50-401N	adult, natural
50-401B	adult, blue
50-402N	pediatric, natural
50-402B	pediatric, blue
50-402G	pediatric, green
50-402R	pediatric, orange
50-402P	pediatric, pink
50-402A	pediatric, assorted pack (2 ea. blue, 2 ea. pink, 1 ea. green, 1 ea. orange)

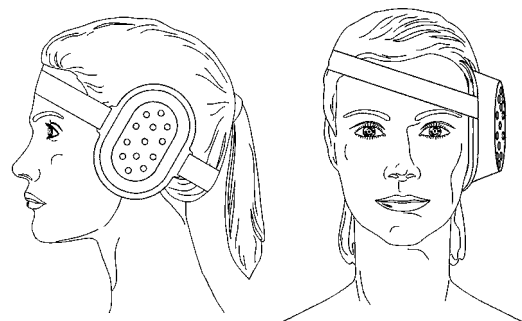



TO PLACE AN ORDER

Our Customer Representatives are available from 8:30 a.m. to 5:00 p.m. (EST) Monday through Friday. Orders can also be placed through our 24-hour fax line, 508-898-2373.

- Customer Service: 800-433-2674 (toll-free in the U.S. and Canada)
- Main Telephone: 508-898-9300
- 24-hour-Fax: 508-898-2373

SHIPMENT

- Standard shipping within the U.S. is FedEx two-day delivery.
- Express orders are accepted until 3:00 p.m. (EST) for next business day delivery nationwide.
- Most orders are shipped the same day.
- Saturday delivery is available in many areas.

WARRANTY

Boston Medical Products, Inc. warrants that reasonable care has been used in the manufacturing of all products. All products are warranted to be free of defects at the time of shipment. For specific product warranties, please refer to individual product instructions for use.

RETURN POLICY

Boston Medical Products, Inc. will accept returned items that are unopened in the original packaging.

- Please call Customer Service to obtain a Return Authorization.
- Returns accepted within 30 days from invoice date.
- A Return Authorization is required for all returns.
- All returns are subject to a restocking fee.
- Merchandise credit only.

